

Power Lifters

Power Lifter 2259M3

Single Column System

Features:

Safe air-operated system featuring an internal, rechargeable 12V battery

Easy to operate:

- Depress the "UP" foot pedal to raise 5,000 cc per hook (total 40 lbs) to a height of 9 feet
- Depress the "DOWN" foot pedal to smoothly lower solution containers and make changeovers
- Heavy cast aluminum base provides stability at maximum and minimum capacity
- 4-hook independent hanger assembly safely holds a total of four 5,000 cc irrigation containers
- Battery operation eliminates hazardous electric cords on O.R. floor

Specifications for Power Lifter

Capacity - Up to four 5000 cc containers (one bag per hook)

Material - Anodized aluminum poles / Stainless steel hanger assembly

Cast aluminum base 6 non-conductive casters

Lift Mechanism - Pneumatic

Height - Min. 5' 7" to Max. 9' from floor to top hanger

Shipping Weight - Base - 48 lbs / Pole Assembly - 8 lbs

Operating Voltage - 12V DC (nom.)

Recharge Power - 115V 60Hz / approx. 180MA




Easy on and off hooks make solution change easy and quick



Locking collar is designed for strength to secure the pole in its raised position. An inner sleeve protects the inner tube against marring



The base features a reservoir to collect fluid and reduce spillage on the O.R. floor

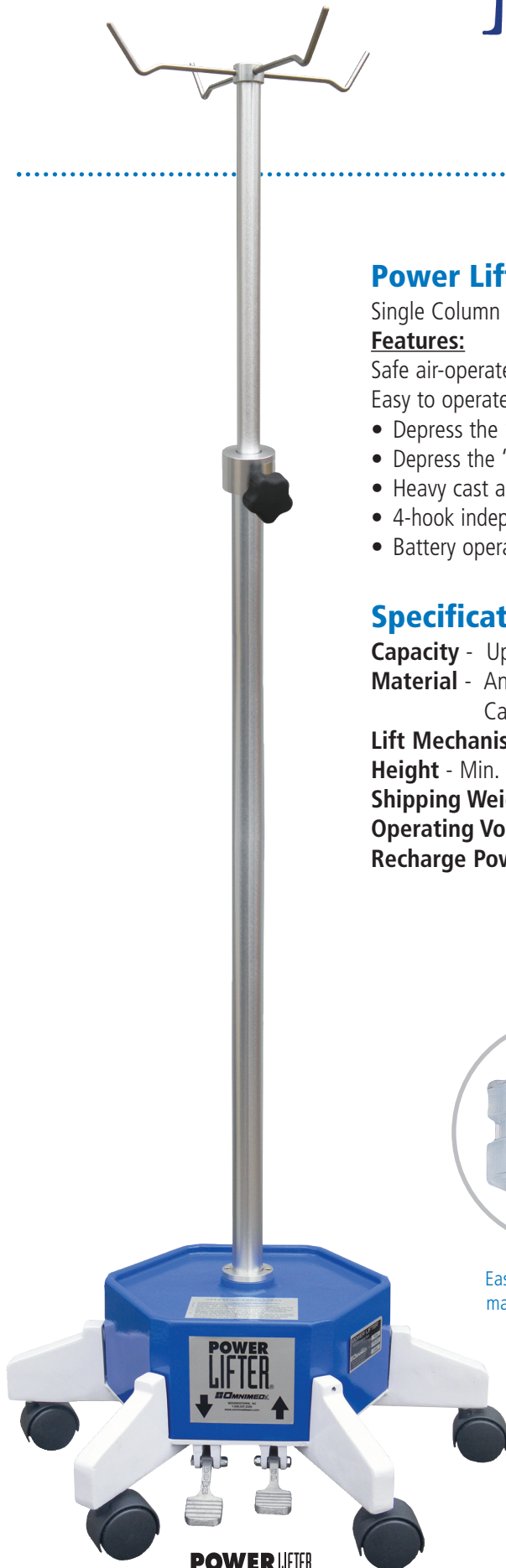

POWER LIFTER.
2259M3