


Factory Installed Options

All options are factory installed for proper fit and finish

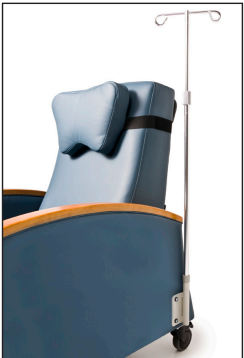

Optional IV Pole Mount is installed on left side of recliner and allows for use of IV Pole.

IV Pole is sold separately.



Optional Side Table is mounted to left side of recliner and can be folded down when not being used.

Specification & Dimensions	Inches	Metric
• Overall Height:	45.5"	(115.6 CM)
• Overall Width:	33"	(83.8 CM)
• Overall Depth:	36.5"	(92.7 CM)
• Width Between Arms:	24"	(61 CM)
• Arm Height From Floor:	26.75"	(67.95 CM)
• Arms Height Above Seat:		
Front of seat	6"	(15.2 CM)
Middle of seat	7.75"	(19.7 CM)
Rear of seat	8.5"	(21.6 CM)
• Seat Height:	21"	(53.3 CM)
• Seat Width:	23"	(58.4 CM)
• Seat Depth:	20.5"	(52.1 CM)
• Depth Fully Reclined:		
1st reclined	64"	(162.6 CM)
Trendelenburg	74"	(188 CM)
• Back Angle Upright:	120°	
Back Angle: Legrest Extended:	130°	
Back Angle: Fully Reclined:	180°	
• Back Height Above Seat	31"	(78.4 CM)
• Recliner Weight:	200 lbs	(91 kg)
• Shipping Weight:	240 lbs	(109 kg)
• Maximum Weight Capacity Evenly Distributed:	350 lbs	(159 kg)

Testing

The Lumex Ortho-Biotic II Recliner meets the following test parameters:

- California Technical Bulletin CA TB133 Flammability Test
- ANSI/BIFMA X5.4-2013 Sections 11,12,14, & 18

Ortho-Biotic II Recliner Warranty

- Frame: Three Years
- Mechanical Components: Three Years
- Upholstered Components: One Year
- Tente Casters: One Year

Subject to Warranty Statement Limitations



Upright Position

Upright Position with Legrest Extended

First Recline Position

Full Recline Position

Trendelenburg Postion

### Ortho-Biotic™ II Recliner

The Ortho-Biotic II Recliner has been made specifically for use in a wide variety of healthcare environments and comes in a beautifully designed and comfortable package. The Ortho-Biotic II Recliner is *proudly Made in the USA!*

The Ortho-Biotic II makes the care, treatment and recovery of patients more comfortable and offers caregivers greater convenience within a wide array of clinical settings.

#### Ideal For:

- Patient Room
- Same Day Surgery
- Post-Anesthesia Recovery
- Oncology/IV Therapy
- Dialysis
- Patient Transport



Ortho-Biotic™ II Recliner with fully upholstered arms



Ortho-Biotic™ II Recliner with wood trimmed arms


### Patient and Caregiver Features

#### Unitized Recliner Mechanism

The Ortho-Biotic II Recliner features a unitized recliner mechanism that allows for independent adjustment of the legrest and back. Easy access Trendelenburg Pedal on right side of the recliner allows for activation of the Trendelenburg Position.



Trendelenburg Pedal

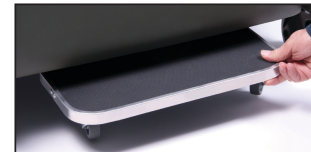

#### Pull Out Patient Foot Plate

For patient comfort and enhanced safety during patient transport, a pull-out footplate is standard. The footplate is easily extended when needed and can be safely stored out of the way when not in use.

#### Legrest-Recline Lever

Conveniently located on the outside of the recliner, the Legrest-Recliner Lever enables the patient or caregiver to easily operate the legrest and adjust the back angle of the recliner.

The Legrest-Recline Lever allows the patient or caregiver to adjust the angle of the back in any position from upright to full recline.



#### Universal Headrest

For enhanced patient comfort in all recline positions, the Ortho-Biotic II comes standard with a Universal Headrest that can easily be adjusted up and down on the back of the recliner. Its unique shape allows the pillow to be used behind the patient's head, neck or back.

#### Composite Foam Seat and Back

For enhanced patient comfort, the seat and back are constructed of a high density foam base with a 1" top layer of visco-elastic foam.



#### Tente® Casters

The Ortho-Biotic II Recliner features four Tente® Swivel Casters specifically designed for healthcare applications. There are three total lock casters with a fourth caster that features a green directional locking tab for improved tracking during patient transport.

#### 46 Fabrics and Colors

- The Lumex Ortho-Biotic II Recliner is offered in a wide array of fabrics and colors. All fabric offerings are specifically designed for healthcare applications and are easy to clean and maintain.
- The Lumex Ortho-Biotic II Recliner meets California Technical Bulletin 133 (TB 133/CA 133) Flammability Test Procedure for Seating Furniture for Use in Public Occupancies
- Please find the most current Lumex Healthcare Seating Swatch Card at [www.grahamfield.com](http://www.grahamfield.com).

#### Standard and Enhanced Fabric Collection

- The standard collection features healthcare grade vinyls that are durable, stain resistant and are available in a wide range of colors to meet any interior design requirements.
- The Enhanced Collection features SILVERGUARD® Silver Ion Technology. SILVERGUARD is a natural environmentally friendly bacteria-fighting defense to inhibit odor and stain causing bacteria, plus the SILVERGUARD technology is EPA registered and FDA approved.